

## Tecnologia

Un compressore Fini MedicAir unisce la semplicità d'uso alla sicurezza di erogare aria compressa priva di impurità.


L'esclusivo essiccatore ad adsorbimento a due colonne garantisce, in ogni condizione di utilizzo, aria secca priva di umidità. Il ciclo di lavoro dell'essiccatore è di circa 5 min., durante i quali una delle due colonne è in fase attiva, trattiene quindi il vapore acqueo trasportato dall'aria compressa, mentre la seconda è in fase di rigenerazione, scaricando tutta l'umidità raccolta attraverso l'elettrovalvola.

Questo ciclo viene invertito tramite un temporizzatore.

## Technology

*A Fini MedicAir compressor is user-friendly and ensures output of pure compressed air, free of impurities.*

*The exclusive twin-tower adsorption drier ensures delivery of dry air, free of moisture, under any conditions of use. The drier work cycle lasts approximately for 5 min.; during this cycle one of the two columns is active and thus holds water vapour conveyed by compressed air, while the second column is in regeneration mode and drains collected moisture through solenoid valve. This cycle is reversed by a timer.*

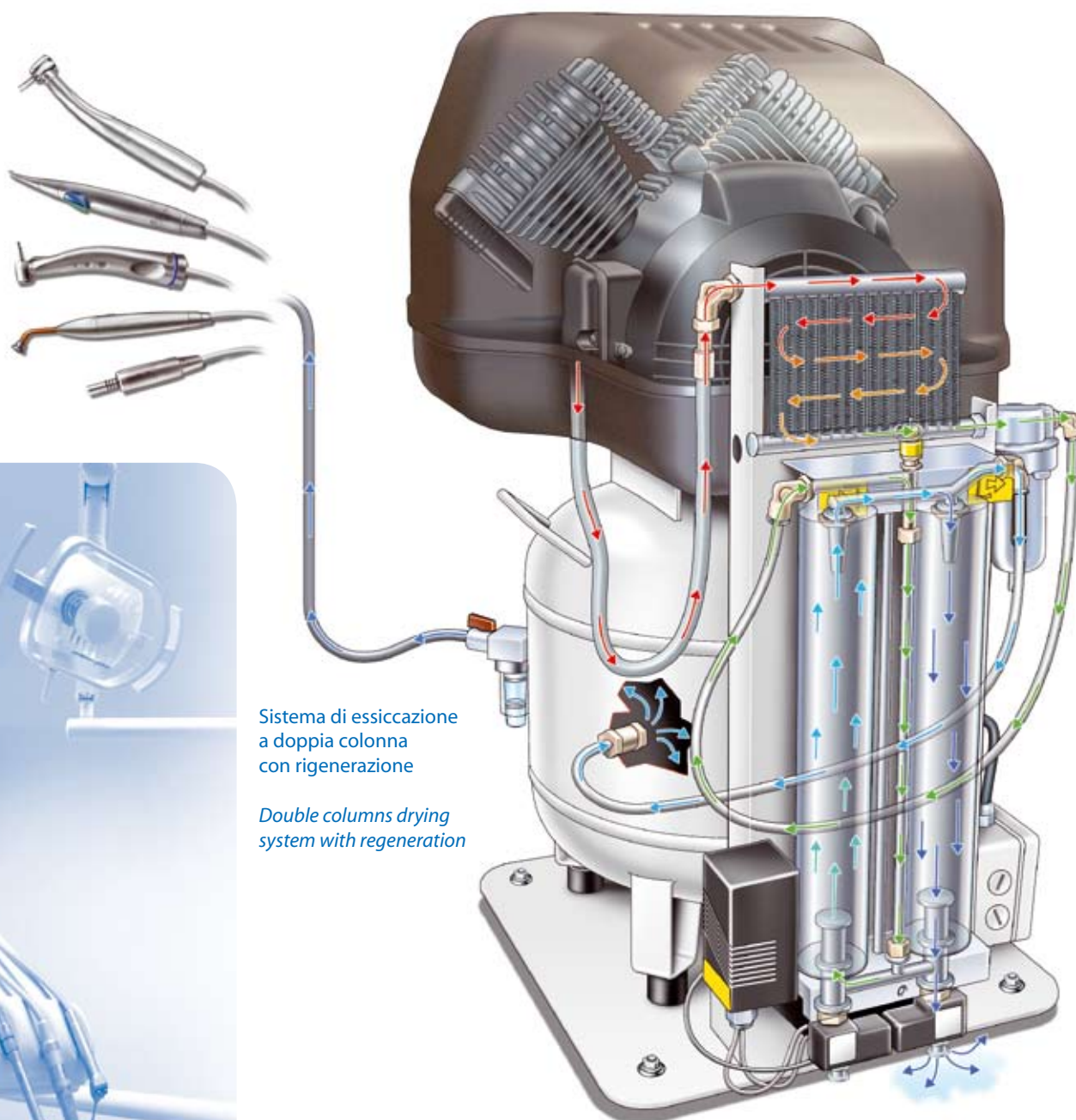

Sistema di essiccazione a doppia colonna con rigenerazione

Double columns drying system with regeneration